



Stent design

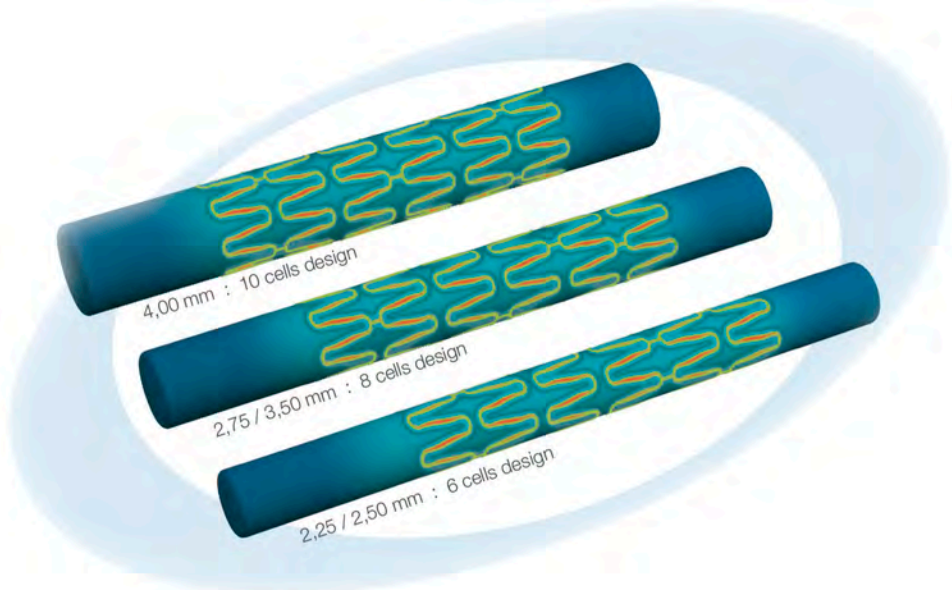

2.25/2.50 - 6 cells design
 2.75/3.50 - 8 cells design
 4.00 and above - 10 cells design
 Homogenous surface distribution.

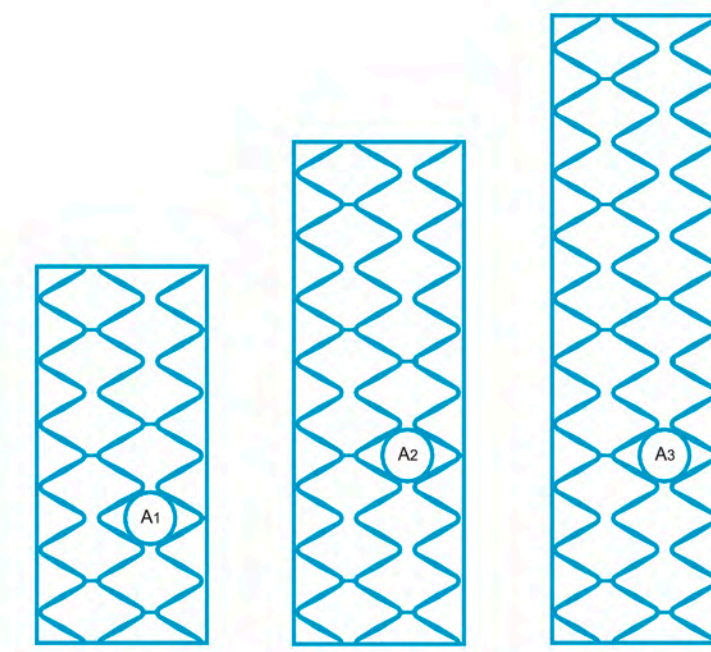

Ordering information

Catalog number	Ø (mm)	L (mm)	Catalog number	Ø (mm)	L (mm)
STCC0001	2.25	8	STCC0025	2.25	24
STCC0002	2.50	8	STCC0026	2.50	24
STCC0003	2.75	8	STCC0027	2.75	24
STCC0004	3.00	8	STCC0028	3.00	24
STCC0005	3.50	8	STCC0029	3.50	24
STCC0006	4.00	8	STCC0030	4.00	24
STCC0031	4.50	8	STCC0039	4.50	24
STCC0032	5.00	8	STCC0040	5.00	24
STCC0007	2.25	12	STCC0041	2.25	28
STCC0008	2.50	12	STCC0042	2.50	28
STCC0009	2.75	12	STCC0043	2.75	28
STCC0010	3.00	12	STCC0044	3.00	28
STCC0011	3.50	12	STCC0045	3.50	28
STCC0012	4.00	12	STCC0046	4.00	28
STCC0033	4.50	12	STCC0047	4.50	28
STCC0034	5.00	12	STCC0048	5.00	28
STCC0013	2.25	16	STCC0049	2.25	32
STCC0014	2.50	16	STCC0050	2.50	32
STCC0015	2.75	16	STCC0051	2.75	32
STCC0016	3.00	16	STCC0052	3.00	32
STCC0017	3.50	16	STCC0053	3.50	32
STCC0018	4.00	16	STCC0054	4.00	32
STCC0035	4.50	16	STCC0055	4.50	32
STCC0036	5.00	16	STCC0056	5.00	32
STCC0019	2.25	20	STCC0057	2.25	36
STCC0020	2.50	20	STCC0058	2.50	36
STCC0021	2.75	20	STCC0059	2.75	36
STCC0022	3.00	20	STCC0060	3.00	36
STCC0023	3.50	20	STCC0061	3.50	36
STCC0024	4.00	20	STCC0062	4.00	36
STCC0037	4.50	20	STCC0063	4.50	36
STCC0038	5.00	20	STCC0064	5.00	36
			STCC0065	2.25	40
			STCC0066	2.50	40
			STCC0067	2.75	40
			STCC0068	3.00	40
			STCC0069	3.50	40
			STCC0070	4.00	40
			STCC0071	4.50	40
			STCC0072	5.00	40

Performing by design


- CroCo Thin
- CroCo Flex
- CroCo Track
- CroCo Performant

Designed by

■ CroCo Thin

■ CroCo Flex

- Extra thin Cobalt Chromium struts 73µm for uncompromising radial strength and radiopacity.
- Unique design, a combination of intermediate and open cells stent for outstanding flexibility and optimal vessel support.

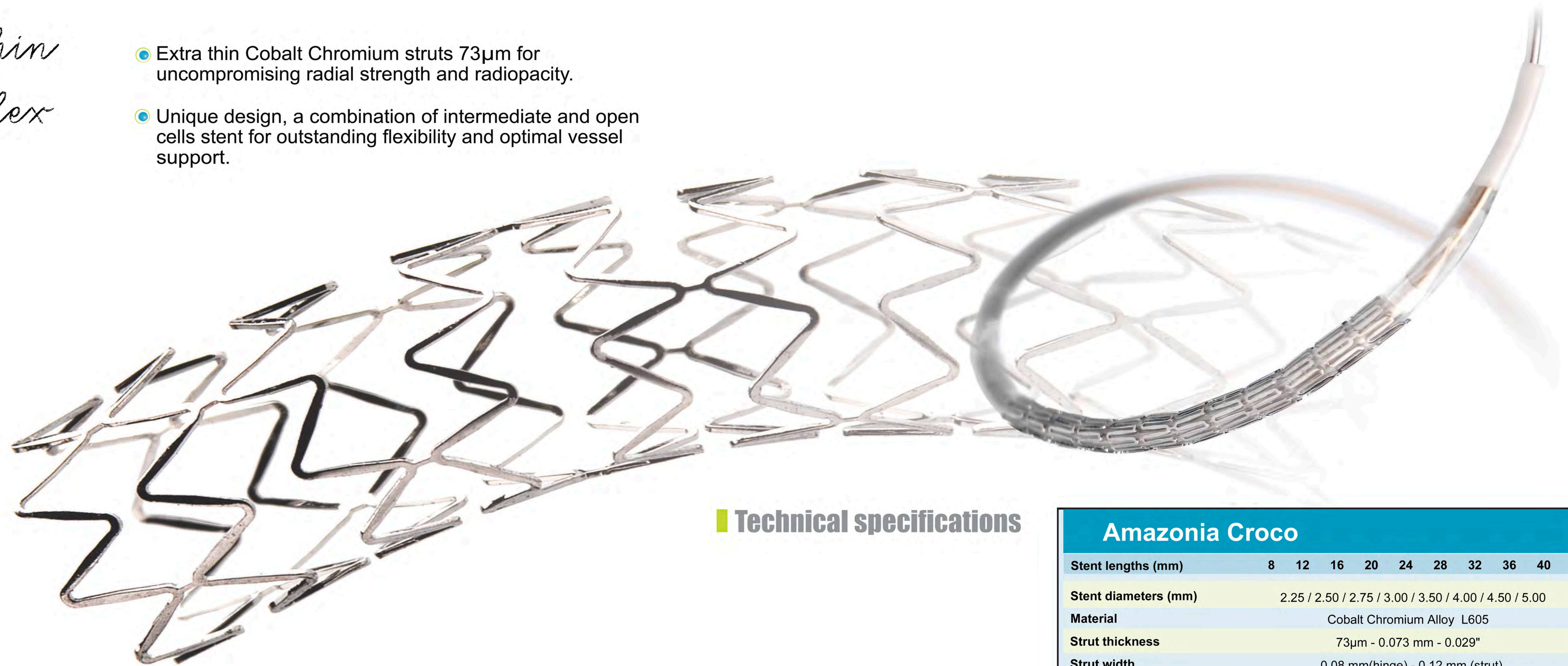

■ Technical specifications

Amazonia CroCo	
Stent lengths (mm)	8 12 16 20 24 28 32 36 40
Stent diameters (mm)	2.25 / 2.50 / 2.75 / 3.00 / 3.50 / 4.00 / 4.50 / 5.00
Material	Cobalt Chromium Alloy L605
Strut thickness	73µm - 0.073 mm - 0.029"
Strut width	0.08 mm(hinge) - 0.12 mm (strut)
Device lengths	Stent length = balloon length = markers distance
Metal to artery ratio	14% average
Nominal pressure	8 ATM
Rated Burst Pressure	14 ATM
Average foreshortening	<1%
Average recoil	<5%
Guiding catheter compatibility	5F (0.056" ID) for all diameters
Guidewire compatibility	0.014" maximum recommended



■ CroCo Track

■ CroCo Performant

- Innovative coated balloon with no overhang and extra low profile for accessing the most tortuous lesions.
- Only designed to perform.